

# ELBOWLIFT™


## SUSPENSION PAD

# Application Guide

STEP

1


Open up the Elbowlift so that the hook-and-loop strap is pulled away from the suspension pad.

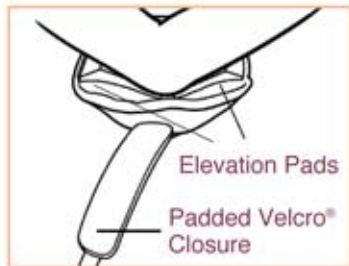

STEP

2


Place tip of patient's elbow over the centering hole with the foam insert pads aligned under the arm and forearm.

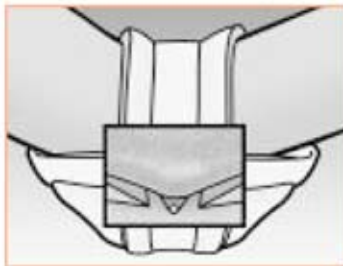

STEP

3

Close the strap across the elbow using the hook and loop closure. Fit the Elbowlift so it is snug but not tight.

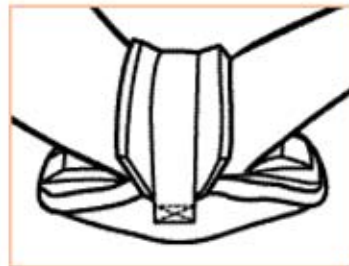

## Care of Elbowlift:

When Elbowlift becomes soiled, hand washing is recommended.

Machine washing Elbowlift in a standard washer is an acceptable alternative to hand washing. The hook-and-loop strap should be closed and placed into a net laundry bag prior to machine washing.

**IMPORTANT:** Remove Elbowlift every eight hours and check the skin or wound status.

For the prevention or treatment of heel ulcers, use



and



Suspension Boots.



Contact details:

Phone: 01986 798120 Fax: 01986 798040

Email: [info@vmorthotics.co.uk](mailto:info@vmorthotics.co.uk) Web: [www.vmorthotics.co.uk](http://www.vmorthotics.co.uk)



© 2001, DM Systems, Inc.  
All rights reserved.